

## Sure Stance™ Knee (THK-4PSC)



### The Only 4-Bar Knee with the Additional Security of Weight-Activated Stance Control.

The Sure Stance Knee's Stumble Control™ feature locks flexion under load from 0 to 35°. Extension remains unlocked allowing the user to smoothly regain balance. Stance Control sensitivity is controlled via a precision adjustment.

#### Primary benefits:

- ✓ **The Sure Stance uniquely combines the proven biomechanical benefits of 4-Bar Swing Phase & the safety net of Weight-Activated Stance Control.**
- ✓ At mid swing-phase, Leg Length is shortened providing increased toe-ground clearance. This drastically reduces both the chance of tripping as well as daily energy consumption.
- ✓ The Sure Stance's Swing Phase Shortening also eliminates the need to shorten the pylon excessively to reduce your Patient's chance of tripping. Correct Leg Length promotes smoother gait and preserves lower back function long-term.
- ✓ High resistance Pneumatic "Hydraulic-Like" dampening provides more precise control of heel rise and extension dampening.

**Suggested K3 L-Codes\*** **L5816** **L5840** **L5845** **L5848** **L5850**

(Frequently used for K2 Patients as well)

**Suggested K2 L-Codes\*** **L5611\*** **GL\*\*** **L5816** **L5845** **L5850**

( See Full Suggested Coding Statement on pg. 53 )

Complete Info & Images @ [daw-usa.com](http://daw-usa.com)

#### Knee Specifications:

Stock Number	THK-4PSC
Functional Level	K3 / K2
Amputation Level	TF or Hip Disarticulation
Weight Limit	250lb (113kg)
Knee Weight	2.4lb (1097g)
Construction Material	Highest Grade Aerospace Alloy
Build Height	(See Next Page)
Stability Controls	4-Bar Geometric Lock plus Weight Activated Stance Control Lock up to 35° flexion
Swing Controls	High Resistance Pneumatic Extension / Flexion Dampening
Maximum Flexion Angle	Approx. 140°
Proximal Connection	Male Pyramid (See Selector Chart, Next Page)
Distal Connection	30mm Tube Clamp
Warranty	2 years included, Upgrade for an Additional 3 Years

#### Order Includes:

- 1x** Practitioner's Manual
- 1x** Suggested L-Codes Letter

#### Recommended Foot:




**K3 Pro-Action™ Foot  
w/ Multi-Axial Ankle**


Full info on page 64

# Sure Stance™ Knee

## Proximal Attachments Selector Chart

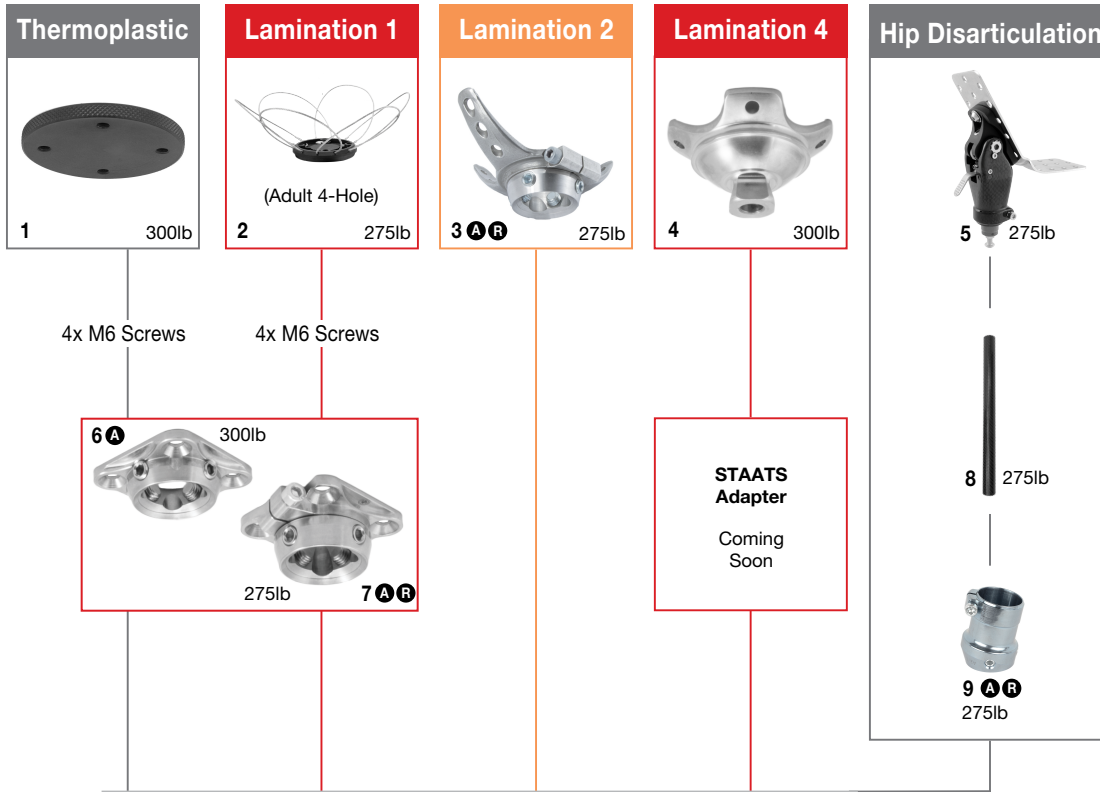

### Build Height

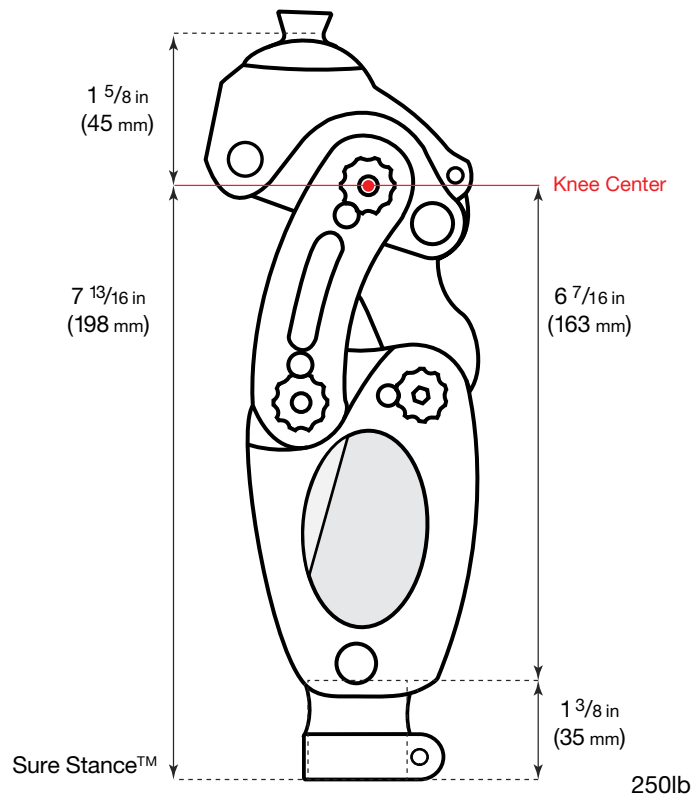

	STOCK	MATERIAL	FULL INFO
1	TSC-T	<b>A</b>	Pg. 75
2	TSC-BX	<b>AA</b>	Pg. 73
3	GUPS-FLAM3PRROT	<b>S</b>	Pg. 74
4	GUPT-MLAM	<b>T</b>	Pg. 73
5	TGH-01	<b>C</b>	Pg. 56
6	GUPT-F4H	<b>T</b>	Pg. 79
7	GUPT-F4HROT	<b>T</b>	Pg. 79
8	TTG-14	<b>C</b>	Pg. 78
9	GUPT-FCLAMP-ANG	<b>T</b>	Pg. 84

MATERIALS
<b>A</b> = Aluminum Alloy
<b>AA</b> = Areospace Aluminum
<b>C</b> = Carbon Graphite
<b>S</b> = Stainless Steel
<b>T</b> = Titanium

ADJUSTABILITY
<b>A</b> = Angular
<b>R</b> = Rotational
<b>S</b> = Sliding