

# SLK™ 3.0 KD Knee

(TGK-5PS10SLKMTX)



PDAC Approved for Microprocessor Swing and Stance!

**Self-Learning Technology:**

The SLK 3.0 KD utilizes Self-Learning Technology to self-program and continually self-adjust as the Amputee’s gait progresses. The Prosthetist simply selects the preferred “Matrix,” then allows the knee to program itself as the Amputee walks.

**The Function of the Microprocessor:**

The SLK 3.0 KD employs 2 types of sensors. One sensor measures the precise flexion/extension speed and angle. The other sensor, a 3-Dimensional gyroscopic sensor, detects the acceleration/deceleration of even the slightest movements in any direction in space. Instantly gathering data from this dual sensor system enables the Microprocessor to accurately recreate natural biomechanical gait based on the unique walking pattern of the Amputee. Contrary to other MPKs, the Self-Learning Knee (SLK) will forever continue adjusting to the progression in the Amputee’s gait.

**Primary benefits:**

- ✓ Provides Cosmetic Evenness in Seated Position.
- ✓ 1-Step Programming.
- ✓ Intelligent Variable Cadence.
- ✓ Natural Polycentric Toe-Clearance in Swing.
- ✓ Custom-ICR™ Stability Adjustment.
- ✓ Ultra-Lite Carbon Fiber Construction.
- ✓ Precise Adjustability of Stance Flexion Stability.



**PDAC Verified L-Codes:** L5613 L5845 L5848 L5856

A copy of this knee’s PDAC Coding Verification Letter can be provided upon request.



**SLK™ 3.0 KD Webpage**


Scan This QR Code For Up-To-Date Info

**Knee Specifications:**

Stock Number	TGK-5PS10SLKMTX
Functional Level	K3 or K4
Amputation Level	Knee Disarticulation
K3 Weight Limit	275lb (125kg)
K4 Weight Limit	220lb (100kg)
Knee Weight	2.7lb (1227g)
Stability Controls	Microprocessor Controlled Stance, Custom-ICR™ Stability Adjustment & Stance Flexion Stability Adjustment
Swing Control	Microprocessor Controlled Swing
Maximum Flexion Angle	Approx. 135°
Proximal Connection	M6 4-Hole Pattern or Single Hole for 7/16 Bolt
Distal Connection	30mm Tube Clamp
Warranty	2 years included, Upgrade for an Additional 3 Years

**Order Includes:**

- 1x Bluetooth Remote
- 1x Charging Cable
- 1x Practitioner’s Manual
- 1x Suggested L-Codes Letter

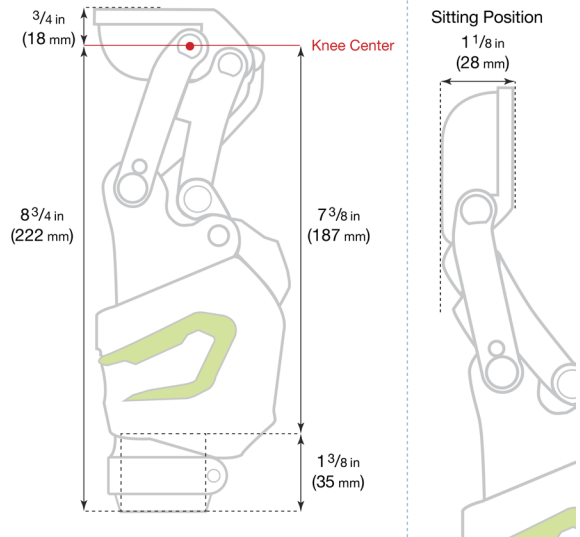

Customer Support (800) 252-2828 • www.daw-usa.com

## Adjustable Features


(Complete Adjustment Instructions Can Be Found In This Knee's Practitioner's Manual, Included With Purchase.)



**Stance Flexion Adjustment**  
(Up To 8°)



**Geometric Stability Adjustment**  
The **Custom-ICR Stability™**  
Adjustment Lets You Dial-In Stability Without Affecting Alignment.

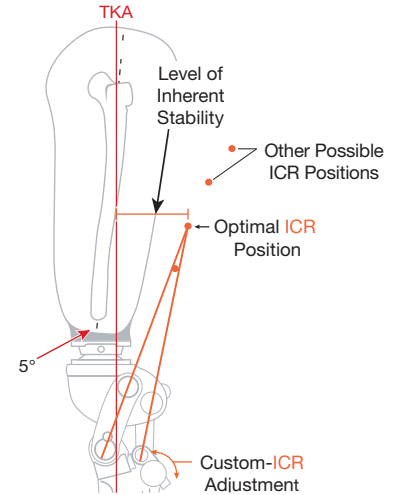

More About Finding Your Patient's Optimal Custom-ICR Stability Setting on page 36.



**1-Step Programming**  
Simply select the preferred "Matrix" using the SLK Remote, then the knee programs itself as the Amputee walks.

## Popular Proximal Options

Proximally Attach Any Adult 4-Hole Connector




**Lo-Pro Rotator™**  
(#: TKR-01)  
Weight Limit: 220lbs



**3-Prong KD Adapter, Stainless Steel**  
(#: TSC-KDL)  
Weight Limit: 275lbs

## Proximal Connector Flow Chart

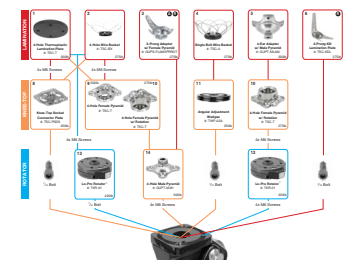

**View All Proximal Options**  
on page 34

## Recommended Foot



**K3 Pro-Action™ Foot w/ Multi-Axial Ankle**  
Weight Limit: up to 300lbs